

Pedicraft Stay Safe Bed

Operations Manual

Before using the Stay Safe Bed, carefully read this Operations Manual.

Be sure you understand all the information in this manual before using the bed.

Instructions for Use	3
Cleaning Instructions	4
Maintenance & Warranty	5
Bed Assembly Instructions	6
Vinyl Top Assembly Instructions	8

For technical assistance call 800-223-7649

Pedicraft, Inc
4134 St. Augustine Road
Jacksonville, FL 32207
800-223-7649
904-396-9627
Fax 904-346-3947
www.pedicraft.com



Orig. - 101414

Note to end user:

Record the following information about your bed and retain this manual for reference.
Please have this information available if you contact Pedicraft about your bed.

Model number _____

Serial number _____

Date received _____

Purchased from _____

Address _____

City, State. Zip _____

Contact name _____

Contact phone _____

Maintenance check (see page 5)

Date: _____	By: _____	Date: _____	By: _____	Date: _____	By: _____
Date: _____	By: _____	Date: _____	By: _____	Date: _____	By: _____
Date: _____	By: _____	Date: _____	By: _____	Date: _____	By: _____
Date: _____	By: _____	Date: _____	By: _____	Date: _____	By: _____
Date: _____	By: _____	Date: _____	By: _____	Date: _____	By: _____
Date: _____	By: _____	Date: _____	By: _____	Date: _____	By: _____
Date: _____	By: _____	Date: _____	By: _____	Date: _____	By: _____
Date: _____	By: _____	Date: _____	By: _____	Date: _____	By: _____

Maintenance Notes

Instructions for Use

The Stay Safe Bed is intended for use during a patient's sleeping hours to prevent falls and wandering while others in the house are sleeping. It is not intended to be used for extended periods when the patient is awake and unsupervised.

Stay Safe Beds are intended for medically fragile patients with chronic illnesses or injuries. They are not intended to restrain a difficult but otherwise healthy child or adult.

- To lower Side Rail: Turn handle and lift slightly, then lower to desired position.
- To raise Side Rail: Hold side rail tube or top rail and lift to desired position. Do not use handle to raise Side Rail.

CAUTION: DO NOT LEAVE PATIENT UNATTENDED WHEN SIDE RAIL IS LOWERED.

- If Bed is equipped with Head Elevation feature: Turn crank at the foot of the bed to raise and lower the head of the bed.

CAUTION: DO NOT LEAVE AN ACTIVE PATIENT UNATTENDED WHEN THE HEAD OF THE BED IS RAISED.

- If Bed is equipped with a Top: Squeeze both knobs and lift rail to open top. Make sure that both latches are secure before releasing the rail.

CAUTION: DO NOT LEAVE PATIENT UNATTENDED WHEN THE TOP IS OPEN.

- The mattress and padding should be routinely inspected for tears in the vinyl cover. Replace the mattress or padding if tears are evident.
- Do not substitute the mattress provided with this bed with any other type or size of mattress. The mattress must fit snugly to prevent gaps which may entrap the patient.

Cleaning of Stay Safe Bed Padding and Sure-Chek[®] Mattress

Wipe clean all parts of the Stay Safe Bed Padding with lukewarm, soapy water on a regular basis.

Clean the Sure-Chek[®] Mattress with mild soap and lukewarm water. Rinse with clean water. Do not use harsh detergents, cleansers or solvents. Daily cleaning of the Mattress is recommended. The mattress should be routinely inspected for tears in the vinyl cover. Replace the mattress if tears are evident.

AIR DRY ONLY. Do not use any form of heat drying. Heat will damage the Padding and Mattress.

Cleaning of Metal Surfaces

Use a damp sponge or soft cloth with a mild soap or detergent solution. Be sure to follow the manufacturer's directions for diluting cleaning solutions. Thoroughly rinse with clean water and dry. Failure to properly rinse and dry will leave a corrosive residue which may cause damage to the metal surfaces or be harmful to people.

Maintenance

The frequency of routine maintenance will vary depending on how active the patient is. A patient who is very active in the bed may cause the bolts to loosen sooner than expected. Record maintenance checks in the space provided on page 2 of this manual.

- All bolts and screws should be checked at least once every six months. Tighten if necessary. **Use caution to not overtighten.**
- Zippers should be routinely inspected to make sure that they are secure.
- The mattress should be inspected for tears in the vinyl cover. Fluid can leak through small tears. The mattress should be replaced as soon as tears are evident.
- The bed padding should be inspected for holes and tears.

Warranty

Pedicraft, Inc. guarantees the Stay Safe Bed frame components to be free of defects in workmanship or materials for a period of three years from the date of shipment to the purchaser. The Stay Safe Bed mattress and bed frame padding and covers are guaranteed to be free of defects in workmanship or materials for a period of ninety (90) days from the date of shipment to the purchaser. In case any purchaser finds a defect in workmanship or materials, promptly notify Pedicraft at 800-223-7649.

The purchaser must return the article as per instructions and Pedicraft will replace, repair or refund the purchase price, at its sole discretion. All replaced or repaired items will be shipped freight prepaid from the factory to the purchaser. Any extra costs resulting from failure to obtain prior shipping instructions from Pedicraft will be borne by the customer.

Any changes or modifications of the equipment by the purchaser without the express written permission of Pedicraft, Inc., will void this warranty. Failure to follow cleaning instructions may void warranty.



Figure 1

Assembly of the bed is easier with 2 people. We recommend clean, rubber gloves be used during assembly to keep the vinyl padding as clean as possible.

Position the Mattress Deck on its side.

Place one End Frame on its side at one end of the Mattress Deck. The padding identifies the top of the frame.



Figure 2

Place the brackets of the Mattress Deck through the End Frame one at a time past the point where the holes line up.

Bring the bottom bracket back through the End Frame and align the holes in the bracket with the holes in the End Frame.



Figure 3

Pull the square tube of the End Frame over the Mattress Deck bracket and align the holes. The fit will be very snug.

Insert four 2" bolts and washers from the outside of the frame and attach nuts. Do not fully tighten.

Repeat to attach the second End Frame.



Figure 4

Attach the Bottom Frame to both End Frames with four 2" bolts and washers.

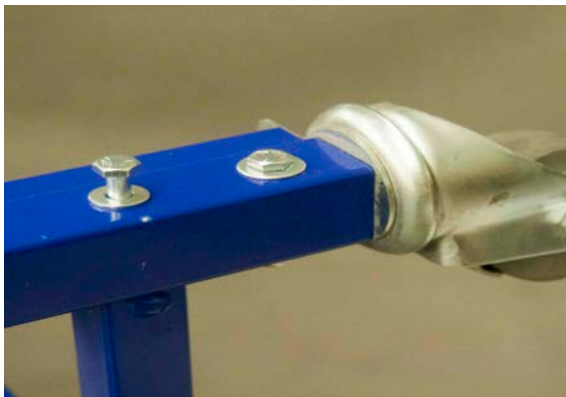


Figure 5

Do not fully tighten.



Figure 6

Stand the Bed in the upright position.

Tighten all bolts starting with the upper right, then lower left and continue diagonally until all bolts are tight (8 Mattress Deck bolts and 4 Bottom Frame bolts).

Use caution to not overtighten.



Figure 7

Remove the nuts from the 3 bolts at the center of each side of the Mattress Deck.

Push the bolts in until they are flush with the Mattress Deck, but do not remove the bolts. Discard the nuts.



Figure 8

Insert one Side Rail into track on each End Frame.

Hold or prop in fully raised position.

Locate Lock Housing at the bottom of the Side Rail.



Figure 9

Grasp Lock Housing and align the 3 holes with the bolts on the Mattress Deck.

Push each of the 3 bolts through the Mattress Deck and thread them into the Lock Housing. Tighten all 3 bolts.

Repeat to attach second Side Rail.



Figure 10

If not adding a Top, assembly is complete.



Figure 11

If your bed is equipped with a Top, proceed to assemble the Top and attach it to the Bed.



Figure 12

We suggest spreading supplied sheet over the box that contained the crib. This provides a work table for assembling the Top.

Place box containing Vinyl Top on top of the large box that contained the bed.



Figure 13

Open the box containing the Top according to the instruction label on the box.

The Vinyl Top parts are packed so that the Top can be assembled in as little space as possible.

Remove the protective foam as you proceed through assembly of the Top.



Figure 14

Start by draping both end sections away from the other parts.

Fold the top side assembly away from the other parts.

Remember to remove the protective foam as you proceed.

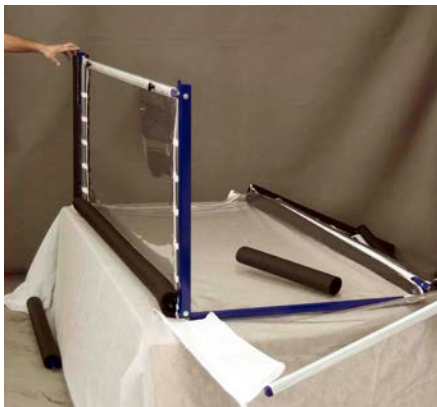


Figure 15

Lift the remaining side assembly upright and align the holes of the side assembly and the blue tube.



Figure 16

Insert 2 1/2" truss head screw with washer through the side assembly into the threaded portion of the blue tube.

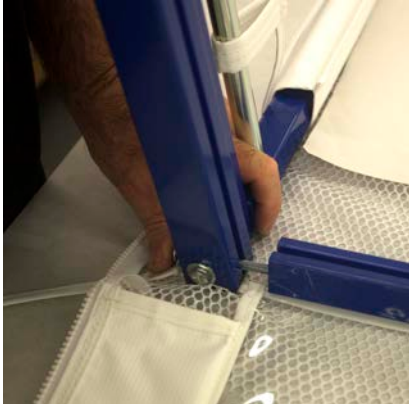


Figure 17



Figure 18

Using the provided no. 3 Phillips screw driver, tighten the screw to bring parts together.

Use caution to not overtighten.



Figure 19

Repeat at the other end of the side assembly.

Use caution to not overtighten.

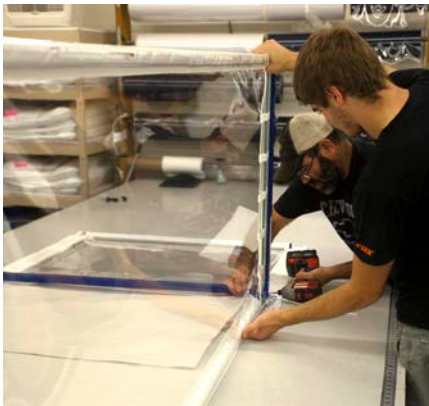


Figure 20

Fold up the other side assembly and attach to the blue tubes in the same fashion.

Use caution to not overtighten.



Figure 21

Raise the end panel on one end of the Top.



Figure 22

As in the previous operation, insert 2½" truss head screw with washer through the side assembly into the threaded portion of the tube of the end panel.

Using the provided no. 3 Phillips screw driver, tighten the screw to bring parts together.

Use caution to not overtighten.



Figure 23

Repeat on the other side of the end panel.

Repeat for the second end panel.

Use caution to not overtighten.



Figure 24

Wrap vinyl sleeve around the tube of the side assembly.



Figure 25

Zip the vinyl sleeve around the tube of the side assembly. Pull the zipper pull completely off of the vinyl sleeve and discard.



Figure 26

Tuck the excess zipper material behind the vinyl sleeve.



Figure 27

Repeat to tuck the excess zipper material behind the vinyl sleeve at the bottom of the Top.

Repeat on the other side of the end panel.

Repeat each step for the opposite end panel.

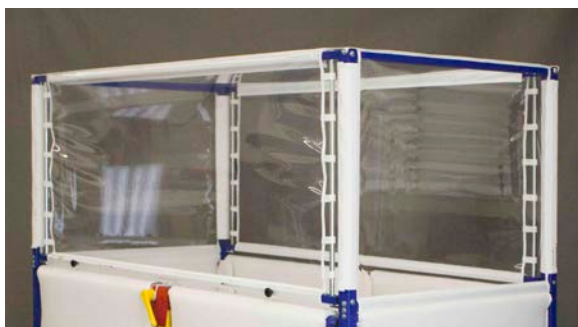


Figure 28

Turn top over (upright).

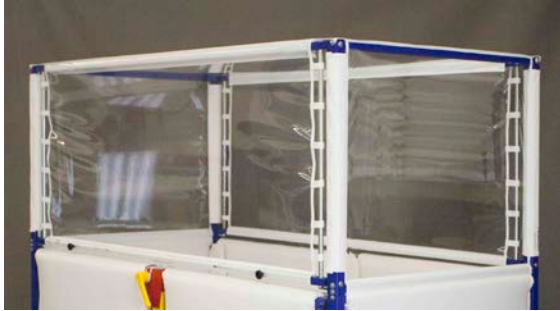


Figure 29

Ensure that black Tube Caps are removed from the corner posts of the Stay Safe Bed.

Lift Top and place on top of the corner posts of the Bed.



Figure 30

Move the Top to align the corner posts of the Top with the corner posts of the Bed.



Figure 31

Turn blue bracket to align holes of bracket with the holes on the Top and Bed.



Figure 32

Insert 1/2" bolt and washer in the top hole of bracket and hand tighten.



Figure 33

Using the provided no. 3 Phillips screw driver, insert 2½" truss head screw (no washer) through the tab of the Top and into the threaded portion of the blue tube of the Bed. Tighten screw.

Use caution to not overtighten.

Tighten both bolts.

Use caution to not overtighten.



Figure 34

Repeat on the other 3 corners of the Top/Bed.

Use caution to not overtighten.

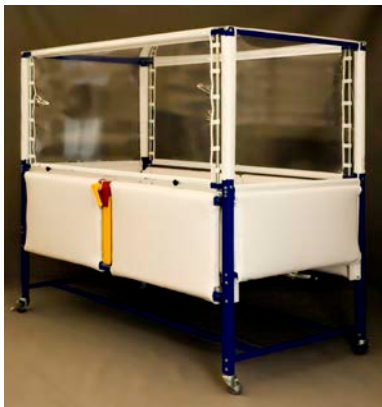


Figure 35

Assembly is complete.